

API Documentation

API Documentation

October 25, 2012

Contents

Contents	1
1 Package pyFAI	2
1.1 Modules	2
1.2 Variables	2
2 Module pyFAI._geometry	3
2.1 Variables	3
3 Module pyFAI.azimuthalIntegrator	4
3.1 Variables	4
3.2 Class AzimuthalIntegrator	4
3.2.1 Methods	5
3.2.2 Properties	20
4 Module pyFAI.bilinear	21
4.1 Variables	21
5 Module pyFAI.detectors	22
5.1 Variables	22
5.2 Class Detector	22
5.2.1 Methods	22
5.2.2 Properties	23
5.3 Class Pilatus	24
5.3.1 Methods	24
5.3.2 Properties	24
5.3.3 Class Variables	25
5.4 Class Pilatus1M	25
5.4.1 Methods	25
5.4.2 Properties	26
5.4.3 Class Variables	26
5.5 Class Pilatus2M	26
5.5.1 Methods	26
5.5.2 Properties	27
5.5.3 Class Variables	27
5.6 Class Pilatus6M	27
5.6.1 Methods	27
5.6.2 Properties	28

5.6.3	Class Variables	28
5.7	Class Fairchild	28
5.7.1	Methods	28
5.7.2	Properties	29
5.8	Class FReLoN	29
5.8.1	Methods	29
5.8.2	Properties	30
6	Module pyFAI.geometry	31
6.1	Variables	31
6.2	Class Geometry	31
6.2.1	Methods	34
6.2.2	Properties	43
7	Module pyFAI.geometryRefinement	44
7.1	Variables	44
7.2	Class GeometryRefinement	44
7.2.1	Methods	44
7.2.2	Properties	47
8	Module pyFAI.histogram	48
8.1	Variables	48
9	Module pyFAI.ocl.azim	49
9.1	Variables	49
10	Module pyFAI.ocl.azim_lut	50
10.1	Variables	50
10.2	Class OCL_LUT_Integrator	50
10.2.1	Methods	50
10.2.2	Properties	51
11	Module pyFAI.peakPicker	52
11.1	Variables	52
11.2	Class PeakPicker	52
11.2.1	Methods	52
11.2.2	Properties	53
11.3	Class ControlPoints	54
11.3.1	Methods	54
11.3.2	Properties	55
11.4	Class Massif	56
11.4.1	Methods	56
11.4.2	Properties	57
12	Module pyFAI.reconstruct	58
12.1	Variables	58
13	Module pyFAI.refinement2D	59
13.1	Variables	59
13.2	Class Refinement2D	59
13.2.1	Methods	59
13.2.2	Properties	60

14 Module pyFAI.relabel	61
14.1 Variables	61
15 Module pyFAI.spline	62
15.1 Variables	62
15.2 Class Spline	62
15.2.1 Methods	62
16 Module pyFAI.splitBBox	67
16.1 Variables	67
17 Module pyFAI.splitBBoxLUT	68
17.1 Variables	68
18 Module pyFAI.splitPixel	69
18.1 Variables	69
19 Module pyFAI.utils	70
19.1 Functions	70
19.2 Variables	72

1 Package pyFAI

1.1 Modules

- **_geometry** (Section 2, p. 3)
- **azimuthalIntegrator** (Section 3, p. 4)
- **bilinear** (Section 4, p. 21)
- **detectors** (Section 5, p. 22)
- **geometry** (Section 6, p. 31)
- **geometryRefinement** (Section 7, p. 44)
- **histogram** (Section 8, p. 48)
- **ocl_azim** (Section 9, p. 49)
- **ocl_azim_lut** (Section 10, p. 50)
- **peakPicker** (Section 11, p. 52)
- **reconstruct** (Section 12, p. 58)
- **refinement2D** (Section 13, p. 59)
- **relabel** (Section 14, p. 61)
- **spline**: This is piece of software aims to manipulate spline files for geometric corrections of the 2D detectors using cubic-spline (Section 15, p. 62)
- **splitBBox** (Section 16, p. 67)
- **splitBBoxLUT** (Section 17, p. 68)
- **splitPixel** (Section 18, p. 69)
- **utils** (Section 19, p. 70)

1.2 Variables

Name	Description
version	Value: '0.7.8'
logger	Value: logging.getLogger("pyFAI.__init__")
__package__	Value: 'pyFAI'

2 Module `pyFAI.geometry`

2.1 Variables

Name	Description
<code>--package--</code>	Value: <code>'pyFAI'</code>
<code>--test--</code>	Value: <code>{}</code>

3 Module *pyFAI.azimuthalIntegrator*

Date: 02/07/2012

Author: Jerome Kieffer

Contact: Jerome.Kieffer@ESRF.eu

Copyright: European Synchrotron Radiation Facility, Grenoble, France

License: GPLv3+

3.1 Variables

Name	Description
<code>__status__</code>	Value: 'beta'
<code>logger</code>	Value: <code>logging.getLogger("pyFAI.azimuthalIntegrator")</code>
<code>ocl</code>	Value: OpenCL devic...
<code>__package__</code>	Value: 'pyFAI'

3.2 Class *AzimuthalIntegrator*



Known Subclasses: `pyFAI.geometryRefinement.GeometryRefinement`

This class is an azimuthal integrator based on P. Boesecke's geometry and histogram algorithm by Manolo S. del Rio and V.A Sole

All geometry calculation are done in the Geometry class

3.2.1 Methods

__init__(*self*, *dist*=1, *poni1*=0, *poni2*=0, *rot1*=0, *rot2*=0, *rot3*=0, *pixel1*=1, *pixel2*=1, *splineFile*=None)

x.__init__(...) initializes *x*; see *help(type(x))* for signature

Parameters

dist: distance sample - detector plan (orthogonal distance, not along the beam), in meter.

poni1: coordinate of the point of normal incidence along the detector's first dimension, in meter

poni2: coordinate of the point of normal incidence along the detector's second dimension, in meter

rot1: first rotation from sample ref to detector's ref, in radians

rot2: second rotation from sample ref to detector's ref, in radians

rot3: third rotation from sample ref to detector's ref, in radians

pixel1: pixel size of the fist dimension of the detector, in meter

pixel2: pixel size of the second dimension of the detector, in meter

splineFile: file containing the geometric distortion of the detector. Overrides the pixel size.

Overrides: *object.__init__*

reset(*self*)

Reset azimuthal integrator in addition to other arrays.

Overrides: *pyFAI.geometry.Geometry.reset*

makeMask(*self*, *data*, *mask*=None, *dummy*=None, *delta_dummy*=None, *invertMask*=None)

Combines a mask

For the mask: 1 for good pixels, 0 for bas pixels

Parameters

data: input array of

mask: input mask

dummy: value of dead pixels

delta_dumy: precision of dummy pixels

invertMask: to force inversion of the input mask

```
xrpd_numpy(self, data, nbPt, filename=None, correctSolidAngle=True, tthRange=None,
mask=None, dummy=None, delta_dummy=None)
```

Calculate the powder diffraction pattern from a set of data, an image. Numpy implementation

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt:	number of points in the output pattern (<i>type=integer</i>)
filename:	file to save data in (<i>type=string</i>)
correctSolidAngle:	if True, the data are divided by the solid angle of each pixel (<i>type=boolean</i>)
tthRange:	The lower and upper range of 2theta. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
mask:	array (same size as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value

Return Value

(2theta, I) in degrees
(*type=2-tuple of 1D arrays*)

```
xrpd_cython(self, data, nbPt, filename=None, correctSolidAngle=True, tthRange=None,
mask=None, dummy=None, delta_dummy=None, pixelSize=None)
```

Calculate the powder diffraction pattern from a set of data, an image. Cython implementation

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt:	number of points in the output pattern (<i>type=integer</i>)
filename:	file to save data in (<i>type=string</i>)
correctSolidAngle:	if True, the data are divided by the solid angle of each pixel (<i>type=boolean</i>)
tthRange:	The lower and upper range of the 2theta. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
mask:	array (same size as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value

Return Value

(2theta, I) in degrees
(*type=2-tuple of 1D arrays*)

xrpd_splitBBox(*self*, *data*, *nbPt*, *filename=None*, *correctSolidAngle=True*, *tthRange=None*, *chiRange=None*, *mask=None*, *dummy=None*, *delta_dummy=None*, *dark=None*, *flat=None*)

Calculate the powder diffraction pattern from a set of data, an image. Cython implementation

TODO: add in the cython part a dark and a flat images to be corrected on the fly. Flat should be combined with solid-angle

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt:	number of points in the output pattern (<i>type=integer</i>)
filename:	file to save data in ascii format 2 column (<i>type=string</i>)
correctSolidAngle:	if True, the data are divided by the solid angle of each pixel (<i>type=boolean</i>)
tthRange:	The lower and upper range of the 2theta. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
chiRange:	The lower and upper range of the chi angle. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
mask:	array (same size as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value
dark:	dark noise image
flat:	flat field image

Return Value

(2theta, I) in degrees
(*type=2-tuple of 1D arrays*)

xrpd_splitPixel(*self*, *data*, *nbPt*, *filename=None*, *correctSolidAngle=True*, *tthRange=None*, *chiRange=None*, *mask=None*, *dummy=None*, *delta_dummy=None*)

Calculate the powder diffraction pattern from a set of data, an image.

Cython implementation

Parameters

data: 2D array from the CCD camera
(*type=ndarray*)

nbPt: number of points in the output pattern
(*type=integer*)

filename: file to save data in
(*type=string*)

mask: array (same size as image) with 0 for masked pixels, and 1 for valid pixels

dummy: value for dead/masked pixels

delta_dummy: precision for dummy value

Return Value

(2theta, I) in degrees
(*type=2-tuple of 1D arrays*)

```
xrpd(self, data, nbPt, filename=None, correctSolidAngle=True, tthRange=None,
chiRange=None, mask=None, dummy=None, delta_dummy=None, dark=None, flat=None)
```

Calculate the powder diffraction pattern from a set of data, an image. Cython implementation

TODO: add in the cython part a dark and a flat images to be corrected on the fly. Flat should be combined with solid-angle

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt:	number of points in the output pattern (<i>type=integer</i>)
filename:	file to save data in ascii format 2 column (<i>type=string</i>)
correctSolidAngle:	if True, the data are divided by the solid angle of each pixel (<i>type=boolean</i>)
tthRange:	The lower and upper range of the 2theta. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
chiRange:	The lower and upper range of the chi angle. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
mask:	array (same size as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value
dark:	dark noise image
flat:	flat field image

Return Value

(2theta, I) in degrees
(*type=2-tuple of 1D arrays*)

```
xrpd_OpenCL(self, data, nbPt, filename=None, correctSolidAngle=True, tthRange=None,
mask=None, dummy=None, delta_dummy=None, devicetype='gpu', useFp64=True,
platformid=None, deviceid=None, safe=True)
```

Calculate the powder diffraction pattern from a set of data, an image. Cython implementation

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt:	number of points in the output pattern (<i>type=integer</i>)
filename:	file to save data in ascii format 2 column (<i>type=string</i>)
correctSolidAngle:	if True, the data are divided by the solid angle of each pixel (<i>type=boolean</i>)
tthRange:	The lower and upper range of the 2theta. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
chiRange:	The lower and upper range of the chi angle. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional, disabled for now</i>)
mask:	array (same size as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value
devicetype:	OpenCL specific parameters: "cpu" or "gpu" or "all" or "def"
useFp64:	shall histogram be done in double precision (advised)
platformid:	platform number
deviceid:	device number
safe:	set to false if you think your GPU is already set-up correctly (2theta, mask, solid angle...)

Return Value

(2theta, I) in degrees
(*type=2-tuple of 1D arrays*)

```
setup_LUT(self, shape, nbPt, mask=None, tthRange=None, chiRange=None)
```

```
xrpd_LUT(self, data, nbPt, filename=None, correctSolidAngle=True, tthRange=None,
chiRange=None, mask=None, dummy=None, delta_dummy=None, safe=True)
```

Calculate the powder diffraction pattern from a set of data, an image. Cython implementation using a Look-Up Table.

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt:	number of points in the output pattern (<i>type=integer</i>)
filename:	file to save data in ascii format 2 column (<i>type=string</i>)
correctSolidAngle:	if True, the data are divided by the solid angle of each pixel (<i>type=boolean</i>)
tthRange:	The lower and upper range of the 2theta. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
chiRange:	The lower and upper range of the chi angle. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional, disabled for now</i>)
mask:	array (same size as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value LUT specific parameters:
safe:	set to false if you think your GPU is already set-up correctly (2theta, mask, solid angle...)

Return Value

(2theta, I) in degrees
(*type=2-tuple of 1D arrays*)

```
xrpd_LUT_OCL(self, data, nbPt, filename=None, correctSolidAngle=True,
tthRange=None, chiRange=None, mask=None, dummy=None, delta_dummy=None,
safe=True, devicetype='all', platformid=None, deviceid=None)
```

Calculate the powder diffraction pattern from a set of data, an image. Cython implementation using a Look-Up Table.

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt:	number of points in the output pattern (<i>type=integer</i>)
filename:	file to save data in ascii format 2 column (<i>type=string</i>)
correctSolidAngle:	if True, the data are divided by the solid angle of each pixel (<i>type=boolean</i>)
tthRange:	The lower and upper range of the 2theta. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
chiRange:	The lower and upper range of the chi angle. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional, disabled for now</i>)
mask:	array (same size as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value LUT specific parameters:
safe:	set to false if you think your GPU is already set-up correctly (2theta, mask, solid angle...) OpenCL specific parameters:
devicetype:	can be "all", "cpu" or "gpu"

Return Value

(2theta, I) in degrees
(*type=2-tuple of 1D arrays*)

```
xrpd2_numpy(self, data, nbPt2Th, nbPtChi=360, filename=None, correctSolidAngle=True,
tthRange=None, chiRange=None, mask=None, dummy=None, delta_dummy=None)
```

Calculate the 2D powder diffraction pattern (2Theta,Chi) from a set of data, an image

Pure numpy implementation (VERY SLOW !!!)

Parameters

data: 2D array from the CCD camera
(*type=ndarray*)

nbPt2Th: number of points in the output pattern in the Radial (horizontal) axis (2 theta)

nbPtChi: number of points in the output pattern along the Azimuthal (vertical) axis (chi)
(*type=integer*)

filename: file to save data in
(*type=string*)

mask: array (same size as image) with 0 for masked pixels, and 1 for valid pixels

dummy: value for dead/masked pixels

delta_dummy: precision for dummy value

nbPt: (*type=integer*)

Return Value

azimuthally regrouped data, 2theta pos and chipos
(*type=3-tuple of ndarrays*)


```
xrpd2_histogram(self, data, nbPt2Th, nbPtChi=360, filename=None,  
correctSolidAngle=True, tthRange=None, chiRange=None, mask=None, dummy=None,  
delta_dummy=None)
```

Calculate the 2D powder diffraction pattern (2Theta,Chi) from a set of data, an image

Cython implementation: fast but incaccurate

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt2Th:	number of points in the output pattern in the Radial (horizontal) axis (2 theta)
nbPtChi:	number of points in the output pattern along the Azimuthal (vertical) axis (chi) (<i>type=integer</i>)
filename:	file to save data in (<i>type=string</i>)
mask:	array (same siza as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value
nbPt:	(<i>type=integer</i>)

Return Value

azimuthaly regrouped data, 2theta pos and chipos
(*type=3-tuple of ndarrays*)

```
xrpd2_splitBBox(self, data, nbPt2Th, nbPtChi=360, filename=None,
correctSolidAngle=True, tthRange=None, chiRange=None, mask=None, dummy=None,
delta_dummy=None)
```

Calculate the 2D powder diffraction pattern (2Theta,Chi) from a set of data, an image

Split pixels according to their coordinate and a bounding box

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt2Th:	number of points in the output pattern in the Radial (horizontal) axis (2 theta)
nbPtChi:	number of points in the output pattern along the Azimuthal (vertical) axis (chi) (<i>type=integer</i>)
filename:	file to save data in (<i>type=string</i>)
correctSolidAngle:	if True, the data are divided by the solid angle of each pixel (<i>type=boolean</i>)
tthRange:	The lower and upper range of the 2theta. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
chiRange:	The lower and upper range of the azimuthal angle. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
mask:	array (same size as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value
nbPt:	(<i>type=integer</i>)

Return Value

azimuthally regrouped data, 2theta pos. and chi pos.
(*type=3-tuple of ndarrays*)

```
xrpd2_splitPixel(self, data, nbPt2Th, nbPtChi=360, filename=None,
correctSolidAngle=True, tthRange=None, chiRange=None, mask=None, dummy=None,
delta_dummy=None)
```

Calculate the 2D powder diffraction pattern (2Theta,Chi) from a set of data, an image

Split pixels according to their corner positions

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt2Th:	number of points in the output pattern in the Radial (horizontal) axis (2 theta)
nbPtChi:	number of points in the output pattern along the Azimuthal (vertical) axis (chi) (<i>type=integer</i>)
filename:	file to save data in (<i>type=string</i>)
correctSolidAngle:	if True, the data are divided by the solid angle of each pixel (<i>type=boolean</i>)
tthRange:	The lower and upper range of the 2theta. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
chiRange:	The lower and upper range of the azimuthal angle. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
mask:	array (same size as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value
nbPt:	(<i>type=integer</i>)

Return Value

azimuthally regrouped data, 2theta pos. and chi pos.
(*type=3-tuple of ndarrays*)

```
xrpd2(self, data, nbPt2Th, nbPtChi=360, filename=None, correctSolidAngle=True,
tthRange=None, chiRange=None, mask=None, dummy=None, delta_dummy=None)
```

Calculate the 2D powder diffraction pattern (2Theta,Chi) from a set of data, an image

Split pixels according to their coordinate and a bounding box

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt2Th:	number of points in the output pattern in the Radial (horizontal) axis (2 theta)
nbPtChi:	number of points in the output pattern along the Azimuthal (vertical) axis (chi) (<i>type=integer</i>)
filename:	file to save data in (<i>type=string</i>)
correctSolidAngle:	if True, the data are divided by the solid angle of each pixel (<i>type=boolean</i>)
tthRange:	The lower and upper range of the 2theta. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
chiRange:	The lower and upper range of the azimuthal angle. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
mask:	array (same size as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value
nbPt:	(<i>type=integer</i>)

Return Value

azimuthally regrouped data, 2theta pos. and chi pos.
(*type=3-tuple of ndarrays*)

```
saxs(self, data, nbPt, filename=None, correctSolidAngle=True, variance=None,
qRange=None, chiRange=None, mask=None, dummy=None, delta_dummy=None,
method='bbox')
```

Calculate the azimuthal integrated Saxs curve

Multi algorithm implementation (tries to be bullet proof)

Parameters

data:	2D array from the CCD camera (<i>type=ndarray</i>)
nbPt:	number of points in the output pattern (<i>type=integer</i>)
filename:	file to save data to (<i>type=string</i>)
correctSolidAngle:	if True, the data are divided by the solid angle of each pixel (<i>type=boolean</i>)
variance:	array containing the variance of the data (<i>type=ndarray</i>)
qRange:	The lower and upper range of the setter vector q. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
chiRange:	The lower and upper range of the chi angle. If not provided, range is simply (data.min(), data.max()). Values outside the range are ignored. (<i>type=(float, float), optional</i>)
mask:	array (same size as image) with 0 for masked pixels, and 1 for valid pixels
dummy:	value for dead/masked pixels
delta_dummy:	precision for dummy value
method:	can be "numpy", "cython", "BBox" or "splitpixel"

Return Value

azimuthally regrouped data, 2theta pos. and chi pos.
(*type=3-tuple of ndarrays*)

```
makeHeaders(self, hdr='#')
```

Return Value

a string to be used for headers

```
save1D(self, filename, dim1, I, error=None, dim1_unit='2th_deg')
```

```
save2D(self, filename, I, dim1, dim2, dim1_unit='2th')
```

Inherited from pyFAI.geometry.Geometry(Section 6.2)

__repr__(), calcfrom1d(), chi(), chiArray(), chi_corner(), cornerArray(), cornerQArray(), del_chia(), del_dssa(), del_qa(), del_ttha(), delta2Theta(), deltaChi(), deltaQ(), diffSolidAngle(), getFit2D(), getPyFAI(), get_chia(), get_correct_solid_angle_for_spline(), get_dist(), get_dssa(), get_pixel1(), get_pixel2(), get_poni1(), get_poni2(), get_qa(), get_rot1(), get_rot2(), get_rot3(), get_spline(), get_splineFile(), get_ttha(), get_wavelength(), load(), oversampleArray(), polarization(), qArray(), qCornerFunct(), qFunction(), read(), save(), setChiDiscAtPi(), setChiDiscAtZero(), setFit2D(), setOversampling(), setPyFAI(), set_chia(), set_correct_solid_angle_for_spline(), set_dist(), set_dssa(), set_pixel1(), set_pixel2(), set_poni1(), set_poni2(), set_qa(), set_rot1(), set_rot2(), set_rot3(), set_spline(), set_splineFile(), set_ttha(), set_wavelength(), sload(), solidAngleArray(), tth(), tth_corner(), twoThetaArray(), write()

Inherited from object

__delattr__(), __format__(), __getattr__(), __hash__(), __new__(), __reduce__(), __reduce_ex__(), __setattr__(), __sizeof__(), __str__(), __subclasshook__()

3.2.2 Properties

Name	Description
<i>Inherited from pyFAI.geometry.Geometry (Section 6.2)</i>	
chia, correct_SA_spline, dist, dssa, pixel1, pixel2, poni1, poni2, qa, rot1, rot2, rot3, spline, splineFile, ttha, wavelength	
<i>Inherited from object</i>	
__class__	

4 Module `pyFAI.bilinear`

Date: 21/12/2011

Author: Jerome Kieffer

Contact: jerome.kieffer@esrf.fr

Copyright: 2011, ESRF

License: GPLv3

4.1 Variables

Name	Description
<code>--package--</code>	Value: 'pyFAI'
<code>--test--</code>	Value: {}

5 Module *pyFAI.detectors*

Date: 12/04/2012

Author: J\xc3\xa9r\xc3\xebme Kieffer

Contact: Jerome.Kieffer@ESRF.eu

Copyright: European Synchrotron Radiation Facility, Grenoble, France

License: GPLv3+

5.1 Variables

Name	Description
<code>__status__</code>	Value: 'beta'
<code>logger</code>	Value: <code>logging.getLogger("pyFAI.detectors")</code>
<code>__package__</code>	Value: 'pyFAI'

5.2 Class Detector

object └─ **pyFAI.detectors.Detector**

Known Subclasses: *pyFAI.detectors.FReLoN*, *pyFAI.detectors.Fairchild*, *pyFAI.detectors.Pilatus*

Generic class representing a 2D detector

5.2.1 Methods

```
__init__(self, pixel1=None, pixel2=None, splineFile=None)
```

`x.__init__(...)` initializes `x`; see `help(type(x))` for signature

Overrides: `object.__init__` `exitit`(inherited documentation)

```
__repr__(self)
```

`repr(x)`

Overrides: `object.__repr__` `exitit`(inherited documentation)


```
get_splineFile(self)
```

```
set_splineFile(self, splineFile)
```

```
get_binning(self)
```

```
set_binning(self, bin_size=(1, 1))
```

```
getPyFAI(self)
```

```
getFit2D(self)
```

```
setPyFAI(self, **kwarg)
```

```
setFit2D(self)
```

```
calc_cartesian_positions(self, d1=None, d2=None)
```

Calculate the position of each pixel center in cartesian coordinate and in meter of a couple of coordinates. The half pixel offset is taken into account here !!!

Parameters

d1: ndarray of dimension 1 or 2 containing the Y pixel positions

d2: ndarray of dimension 1 or 2 containing the X pixel positions

Return Value

2-arrays of same shape as d1 & d2 with the position in meter

d1 and d2 must have the same shape, returned array will have the same shape.

```
get_mask(self)
```

Should return a generic mask for the detector

Inherited from object

```
__delattr__(), __format__(), __getattr__(), __hash__(), __new__(), __reduce__(), __reduce_ex__(),
__setattr__(), __sizeof__(), __str__(), __subclasshook__()
```


5.2.2 Properties

Name	Description
splineFile	

continued on next page

Name	Description
binning	
<i>Inherited from object</i>	
__class__	

5.3 Class Pilatus

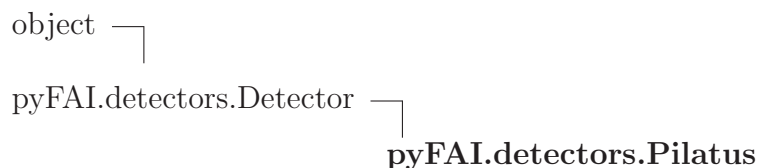

Known Subclasses: pyFAI.detectors.Pilatus1M, pyFAI.detectors.Pilatus2M, pyFAI.detectors.Pilatus6M

Pilatus detector: generic description

5.3.1 Methods

```
__init__(self, pixel1=0.000172, pixel2=0.000172)
x.__init__(...) initializes x; see help(type(x)) for signature
Overrides: object.__init__ extit(inherited documentation)
```

```
get_mask(self)
Returns a generic mask for Pilatus detecors...
Overrides: pyFAI.detectors.Detector.get_mask
```

Inherited from pyFAI.detectors.Detector(Section 5.2)

```
__repr__(), calc_cartesian_positions(), getFit2D(), getPyFAI(), get_binning(), get_splineFile(),
setFit2D(), setPyFAI(), set_binning(), set_splineFile()
```

Inherited from object

```
__delattr__(), __format__(), __getattr__(), __hash__(), __new__(), __reduce__(), __reduce_ex__(),
__setattr__(), __sizeof__(), __str__(), __subclasshook__()
```

5.3.2 Properties

Name	Description
<i>Inherited from pyFAI.detectors.Detector (Section 5.2)</i>	


continued on next page

Name	Description
binning, splineFile	
<i>Inherited from object</i>	
__class__	

5.3.3 Class Variables

Name	Description
MODULE_SIZE	Value: (195, 487)
MODULE_GAP	Value: (17, 7)

5.4 Class *Pilatus1M*

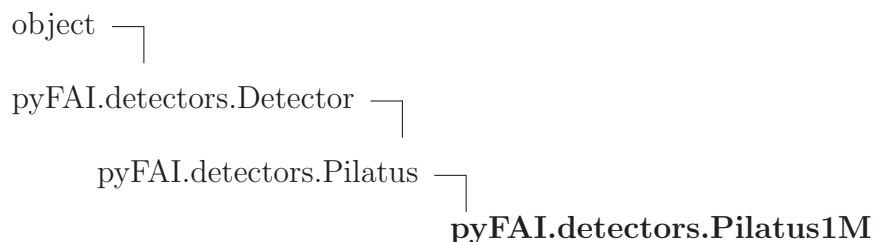

Pilatus 1M detector

5.4.1 Methods

__init__ (<i>self</i> , <i>pixel1</i> =0.000172, <i>pixel2</i> =0.000172) <i>x</i> . __init__ (...) initializes <i>x</i> ; see help(type(<i>x</i>)) for signature Overrides: object. __init__ extit(inherited documentation)
--

Inherited from pyFAI.detectors.Pilatus(Section 5.3)

get_mask()

Inherited from pyFAI.detectors.Detector(Section 5.2)

__repr__(), **calc_cartesian_positions**(), **getFit2D**(), **getPyFAI**(), **get_binning**(), **get_splineFile**(), **setFit2D**(), **setPyFAI**(), **set_binning**(), **set_splineFile**()

Inherited from object

__delattr__(), **__format__**(), **__getattr__**(), **__hash__**(), **__new__**(), **__reduce__**(), **__reduce_ex__**(), **__setattr__**(), **__sizeof__**(), **__str__**(), **__subclasshook__**()


5.4.2 Properties

Name	Description
	<i>Inherited from <code>pyFAI.detectors.Detector</code> (Section 5.2)</i>
	binning, splineFile
	<i>Inherited from object</i>
	<code>__class__</code>

5.4.3 Class Variables

Name	Description
	<i>Inherited from <code>pyFAI.detectors.Pilatus</code> (Section 5.3)</i>
	MODULE_GAP, MODULE_SIZE

5.5 Class *Pilatus2M*

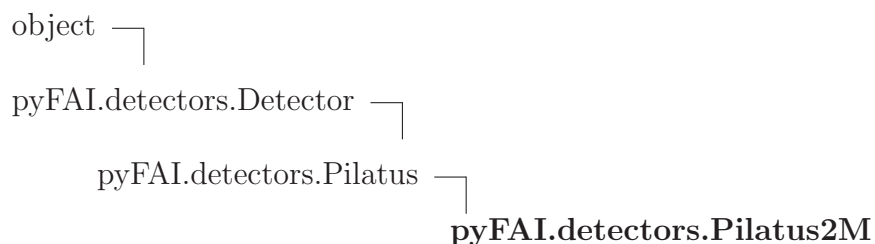

Pilatus 2M detector

5.5.1 Methods

<code>__init__(self, pixel1=0.000172, pixel2=0.000172)</code> <code>x.__init__(...)</code> initializes x; see <code>help(type(x))</code> for signature Overrides: <code>object.__init__</code> <code>extit</code> (inherited documentation)

Inherited from `pyFAI.detectors.Pilatus` (Section 5.3)

`get_mask()`

Inherited from `pyFAI.detectors.Detector` (Section 5.2)

`__repr__()`, `calc_cartesian_positions()`, `getFit2D()`, `getPyFAI()`, `get_binning()`, `get_splineFile()`,
`setFit2D()`, `setPyFAI()`, `set_binning()`, `set_splineFile()`

Inherited from object

`__delattr__()`, `__format__()`, `__getattr__()`, `__hash__()`, `__new__()`, `__reduce__()`, `__reduce_ex__()`,
`__setattr__()`, `__sizeof__()`, `__str__()`, `__subclasshook__()`


5.5.2 Properties

Name	Description
<i>Inherited from <code>pyFAI.detectors.Detector</code> (Section 5.2)</i>	
<code>binning</code> , <code>splineFile</code>	
<i>Inherited from object</i>	
<code>__class__</code>	

5.5.3 Class Variables

Name	Description
<i>Inherited from <code>pyFAI.detectors.Pilatus</code> (Section 5.3)</i>	
<code>MODULE_GAP</code> , <code>MODULE_SIZE</code>	

5.6 Class *Pilatus6M*

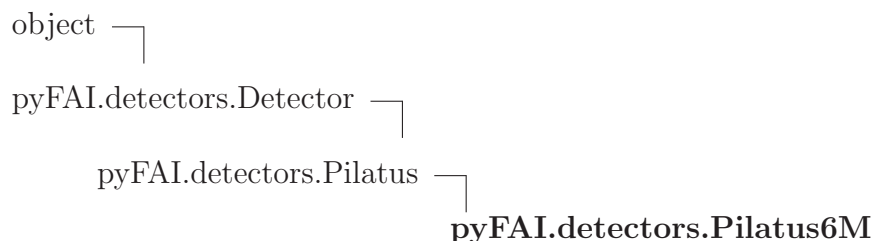

Pilatus 6M detector

5.6.1 Methods

```

__init__(self, pixel1=0.000172, pixel2=0.000172)
x.__init__(...) initializes x; see help(type(x)) for signature
Overrides: object.__init__ extit(inherited documentation)

```

Inherited from `pyFAI.detectors.Pilatus` (Section 5.3)

`get_mask()`

Inherited from `pyFAI.detectors.Detector` (Section 5.2)

`__repr__()`, `calc_cartesian_positions()`, `getFit2D()`, `getPyFAI()`, `get_binning()`, `get_splineFile()`,
`setFit2D()`, `setPyFAI()`, `set_binning()`, `set_splineFile()`

Inherited from object

`__delattr__()`, `__format__()`, `__getattr__()`, `__hash__()`, `__new__()`, `__reduce__()`, `__reduce_ex__()`,
`__setattr__()`, `__sizeof__()`, `__str__()`, `__subclasshook__()`


5.6.2 Properties

Name	Description
<i>Inherited from pyFAI.detectors.Detector (Section 5.2)</i>	
<code>binning</code> , <code>splineFile</code>	
<i>Inherited from object</i>	
<code>__class__</code>	

5.6.3 Class Variables

Name	Description
<i>Inherited from pyFAI.detectors.Pilatus (Section 5.3)</i>	
<code>MODULE_GAP</code> , <code>MODULE_SIZE</code>	

5.7 Class Fairchild

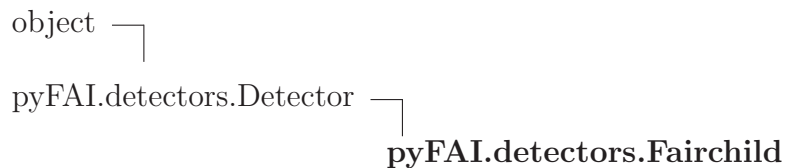

Fairchild Condor 486:90 detector

5.7.1 Methods

`__init__(self, pixel1=1.5e-05, pixel2=1.5e-05)`
`x.__init__(...)` initializes x; see `help(type(x))` for signature
 Overrides: `object.__init__` `extit` (inherited documentation)

Inherited from pyFAI.detectors.Detector (Section 5.2)

`__repr__()`, `calc_cartesian_positions()`, `getFit2D()`, `getPyFAI()`, `get_binning()`, `get_mask()`,
`get_splineFile()`, `setFit2D()`, `setPyFAI()`, `set_binning()`, `set_splineFile()`


Inherited from object

`__delattr__()`, `__format__()`, `__getattr__()`, `__hash__()`, `__new__()`, `__reduce__()`, `__reduce_ex__()`,
`__setattr__()`, `__sizeof__()`, `__str__()`, `__subclasshook__()`

5.7.2 Properties

Name	Description
<i>Inherited from pyFAI.detectors.Detector (Section 5.2)</i>	
<code>binning</code> , <code>splineFile</code>	
<i>Inherited from object</i>	
<code>__class__</code>	

5.8 Class FReLoN

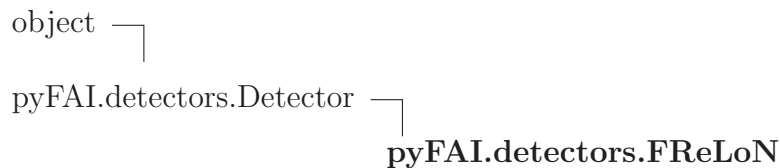

FReLoN detector (spline mandatory to correct for geometric distortion)

5.8.1 Methods

`__init__(self, splineFile)`
`x.__init__(...)` initializes x; see `help(type(x))` for signature
 Overrides: `object.__init__` extit(inherited documentation)

Inherited from pyFAI.detectors.Detector(Section 5.2)

`__repr__()`, `calc_cartesian_positions()`, `getFit2D()`, `getPyFAI()`, `get_binning()`, `get_mask()`,
`get_splineFile()`, `setFit2D()`, `setPyFAI()`, `set_binning()`, `set_splineFile()`

Inherited from object

`__delattr__()`, `__format__()`, `__getattr__()`, `__hash__()`, `__new__()`, `__reduce__()`, `__reduce_ex__()`,
`__setattr__()`, `__sizeof__()`, `__str__()`, `__subclasshook__()`

5.8.2 Properties

Name	Description
	<i>Inherited from pyFAI.detectors.Detector (Section 5.2)</i>
binning, splineFile	
	<i>Inherited from object</i>
__class__	

6 Module *pyFAI.geometry*

Date: 09/06/2012

Author: Jerome Kieffer

Contact: Jerome.Kieffer@ESRF.eu

Copyright: European Synchrotron Radiation Facility, Grenoble, France

License: GPLv3+

6.1 Variables

Name	Description
<code>--status--</code>	Value: 'beta'
<code>logger</code>	Value: <code>logging.getLogger("pyFAI.geometry")</code>
<code>--package--</code>	Value: 'pyFAI'

6.2 Class Geometry

object └─ `pyFAI.geometry.Geometry`

Known Subclasses: `pyFAI.azimuthalIntegrator.AzimuthalIntegrator`

This class is an azimuthal integrator based on P. Boesecke's geometry and histogram algorithm by Manolo S. del Rio and V.A Sole

Detector is assumed to be corrected from "raster orientation" effect.

It is not addressed here but rather in the Detector object or at read time.

Considering there is no tilt:

Detector fast dimension (dim2) is supposed to be horizontal (dimension X of the image)

Detector slow dimension (dim1) is supposed to be vertical, upwards (dimension Y of the image)

The third dimension is chose such as the referential is orthonormal, so dim3 is along image

Demonstration of the equation done using Mathematica.

=====

Axis 1 is along first dimension of detector (when not tilted), this is the slow dimension

```

x1={1,0,0}
Axis 2 is along second dimension of detector (when not tilted), this is the fast dimension
x2={0,1,0}
Axis 3 is along the incident X-Ray beam
x3={0,0,1}
We define the 3 rotation around axis 1, 2 and 3:
rotM1 = RotationMatrix[rot1,x1] = {{1,0,0},{0,cos[rot1],-sin[rot1]},{0,sin[rot1],cos[rot1]}}
rotM2 = RotationMatrix[rot2,x2] = {{cos[rot2],0,sin[rot2]},{0,1,0},{-sin[rot2],0,cos[rot2]}}
rotM3 = RotationMatrix[rot3,x3] = {{cos[rot3],-sin[rot3],0},{sin[rot3],cos[rot3],0},{0,0,1}}

Rotations of the detector are applied first Rot around axis 1, then axis 2 and finally axis 3:
R = rotM3.rotM2.rotM1
R = {{cos[rot2] cos[rot3],cos[rot3] sin[rot1] sin[rot2]-cos[rot1] sin[rot3],cos[rot1] cos[rot3] sin[rot2]+sin[rot1] sin[rot3],-cos[rot2] sin[rot3],cos[rot1] cos[rot3]+sin[rot1] sin[rot2] sin[rot3],-cos[rot1] sin[rot2] sin[rot3],sin[rot2] sin[rot3],cos[rot2] sin[rot1],cos[rot1] cos[rot2]}}
In Python notation:
PForm[R.x1] = [cos(rot2)*cos(rot3),cos(rot2)*sin(rot3),-sin(rot2)]
PForm[R.x2] = [cos(rot3)*sin(rot1)*sin(rot2) - cos(rot1)*sin(rot3),cos(rot1)*cos(rot3) + sin(rot1)*sin(rot2)*sin(rot3),-cos(rot1)*cos(rot2)*sin(rot3)]
PForm[R.x3] = [cos(rot1)*cos(rot3)*sin(rot2) + sin(rot1)*sin(rot3),-(cos(rot3)*sin(rot1)),cos(rot1)*cos(rot2)]

* Coordinates of the Point of Normal Incidence:
PONI = R.{0,0,L}
PForm[PONI] = [L*(cos(rot1)*cos(rot3)*sin(rot2) + sin(rot1)*sin(rot3)),
               L*(-cos(rot3)*sin(rot1)) + cos(rot1)*sin(rot2)*sin(rot3),L*cos(rot1)*cos(rot2)]

* Any pixel on detector plan at coordinate (d1, d2) in meters. Detector is at z=L
P={d1,d2,L}
PForm[R.P] = [t1, t2, t3] =
               [d1*cos(rot2)*cos(rot3) + d2*(cos(rot3)*sin(rot1)*sin(rot2) - cos(rot1)*sin(rot2)*sin(rot3)),
               d1*cos(rot2)*sin(rot3) + d2*(cos(rot1)*cos(rot3) + sin(rot1)*sin(rot2)*sin(rot3)),
               d2*cos(rot2)*sin(rot1) - d1*sin(rot2) + L*cos(rot1)*cos(rot2)]

* Distance sample (origin) to detector point (d1,d2)
FForm[Norm[R.P]] = sqrt(pow(Abs(L*cos(rot1)*cos(rot2) + d2*cos(rot2)*sin(rot1) - d1*sin(rot2)*sin(rot3)),2) +
                        pow(Abs(d1*cos(rot2)*cos(rot3) + d2*(cos(rot3)*sin(rot1)*sin(rot2) - cos(rot1)*sin(rot2)*sin(rot3)),2) +
                        L*(cos(rot1)*cos(rot3)*sin(rot2) + sin(rot1)*sin(rot3)),2) +
                        pow(Abs(d1*cos(rot2)*sin(rot3) + L*(-cos(rot3)*sin(rot1)) + cos(rot1)*sin(rot2)*sin(rot3)),2))

* cos(2theta) is defined as (R.P component along x3) over the distance from origin to detector point
tth = ArcCos [- (R.P).x3/Norm[R.P]]
FForm[tth] = Arccos((-L*cos(rot1)*cos(rot2)) - d2*cos(rot2)*sin(rot1) + d1*sin(rot2))/Norm[R.P]

```

```

sqrt(pow(Abs(L*cos(rot1)*cos(rot2) + d2*cos(rot2)*sin(rot1) - d1*sin
pow(Abs(d1*cos(rot2)*cos(rot3) + d2*(cos(rot3)*sin(rot1)*sin(rot2)
L*(cos(rot1)*cos(rot3)*sin(rot2) + sin(rot1)*sin(rot3))),2) +
pow(Abs(d1*cos(rot2)*sin(rot3) + L*(-(cos(rot3)*sin(rot1)) + cos(r
d2*(cos(rot1)*cos(rot3) + sin(rot1)*sin(rot2)*sin(rot3))),2)))

```

```

* tan(2theta) is defined as sqrt(t1**2 + t2**2) / t3
tth = ArcTan2 [sqrt(t1**2 + t2**2) , t3 ]

```

Getting 2theta from it's tangeant seems both more precise (around beam stop very far from) Currently there is a swich in the method to follow one path or the other.

```

* Tangeant of angle chi is defined as (R.P component along x1) over (R.P component along
chi = ArcTan[((R.P).x1) / ((R.P).x2)]
FForm[chi] = ArcTan2(d1*cos(rot2)*cos(rot3) + d2*(cos(rot3)*sin(rot1)*sin(rot2) - cos(r
L*(cos(rot1)*cos(rot3)*sin(rot2) + sin(rot1)*sin(rot3))),
d1*cos(rot2)*sin(rot3) + L*(-(cos(rot3)*sin(rot1)) + cos(rot1)*sin
d2*(cos(rot1)*cos(rot3) + sin(rot1)*sin(rot2)*sin(rot3)))

```

6.2.1 Methods

```
__init__(self, dist=1, poni1=0, poni2=0, rot1=0, rot2=0, rot3=0, pixel1=1, pixel2=1, splineFile=None, detector=None)
```

x.__init__(...) initializes x; see help(type(x)) for signature

Parameters

dist: distance sample - detector plan (orthogonal distance, not along the beam), in meter.

poni1: coordinate of the point of normal incidence along the detector's first dimension, in meter

poni2: coordinate of the point of normal incidence along the detector's second dimension, in meter

rot1: first rotation from sample ref to detector's ref, in radians

rot2: second rotation from sample ref to detector's ref, in radians

rot3: third rotation from sample ref to detector's ref, in radians

pixel1: pixel size of the first dimension of the detector, in meter

pixel2: pixel size of the second dimension of the detector, in meter

splineFile: file containing the geometric distortion of the detector. Overrides the pixel size.

Overrides: object.__init__

```
__repr__(self)
```

repr(x)

Overrides: object.__repr__ exitit(inherited documentation)

tth(*self*, *d1*, *d2*, *param*=None, *path*='cython')

Calculates the 2theta value for the center of a given pixel (or set of pixels)

Parameters

d1: position(s) in pixel in first dimension (c order)

(type=scalar or array of scalar)

d2: position(s) in pixel in second dimension (c order)

(type=scalar or array of scalar)

path: can be "cos", "tan" or "cython" @return 2theta in radians

Return Value

float or array of floats.

qFunction(*self*, *d1*, *d2*, *param*=None, *path*='cython')

Calculates the q value for the center of a given pixel (or set of pixels) in nm⁻¹

$$q = 4\pi/\lambda \sin(2\theta / 2)$$

Parameters

d1: position(s) in pixel in first dimension (c order)

(type=scalar or array of scalar)

d2: position(s) in pixel in second dimension (c order)

(type=scalar or array of scalar @return q in in nm⁻¹)

Return Value

float or array of floats.

qArray(*self*, *shape*)

Generate an array of the given shape with q(i,j) for all elements.

qCornerFunc(*self*, *d1*, *d2*)

calculate the q_vector for any pixel corner

tth_corner(*self*, *d1*, *d2*)

Calculates the 2theta value for the corner of a given pixel (or set of pixels)

Parameters

d1: position(s) in pixel in first dimension (c order)

(type=scalar or array of scalar)

d2: position(s) in pixel in second dimension (c order)

(type=scalar or array of scalar @return 2theta in radians)

Return Value

float or array of floats.

twoThetaArray(*self*, *shape*)

Generate an array of the given shape with two-theta(i,j) for all elements.

chi(*self*, *d1*, *d2*, *path*='cython')

Calculate the chi (azimuthal angle) for the centre of a pixel at coordinate d1,d2 which in the lab ref has coordinate: $X1 = p1 \cdot \cos(\text{rot2}) \cdot \cos(\text{rot3}) + p2 \cdot (\cos(\text{rot3}) \cdot \sin(\text{rot1}) \cdot \sin(\text{rot2}) - \cos(\text{rot1}) \cdot \sin(\text{rot3})) - L \cdot (\cos(\text{rot1}) \cdot \cos(\text{rot3}) \cdot \sin(\text{rot2}) + \sin(\text{rot1}) \cdot \sin(\text{rot3}))$ $X2 = p1 \cdot \cos(\text{rot2}) \cdot \sin(\text{rot3}) - L \cdot (-\cos(\text{rot3}) \cdot \sin(\text{rot1})) + \cos(\text{rot1}) \cdot \sin(\text{rot2}) \cdot \sin(\text{rot3})) + p2 \cdot (\cos(\text{rot1}) \cdot \cos(\text{rot3}) + \sin(\text{rot1}) \cdot \sin(\text{rot2}) \cdot \sin(\text{rot3}))$ $X3 = -(L \cdot \cos(\text{rot1}) \cdot \cos(\text{rot2})) + p2 \cdot \cos(\text{rot2}) \cdot \sin(\text{rot1}) - p1 \cdot \sin(\text{rot2})$ hence $\tan(\text{Chi}) = X2 / X1$

Parameters

d1: pixel coordinate along the 1st dimension (C convention)

(type=float or array of them)

d2: pixel coordinate along the 2nd dimension (C convention)

(type=float or array of them)

path: can be "tan" (i.e via numpy) or "cython"

Return Value

chi, the azimuthal angle in rad

chi_corner(*self*, *d1*, *d2*)

Calculate the chi (azimuthal angle) for the corner of a pixel at coordinate d1,d2 which in the lab ref has coordinate: $X1 = p1 \cdot \cos(\text{rot2}) \cdot \cos(\text{rot3}) + p2 \cdot (\cos(\text{rot3}) \cdot \sin(\text{rot1}) \cdot \sin(\text{rot2}) - \cos(\text{rot1}) \cdot \sin(\text{rot3})) - L \cdot (\cos(\text{rot1}) \cdot \cos(\text{rot3}) \cdot \sin(\text{rot2}) + \sin(\text{rot1}) \cdot \sin(\text{rot3}))$ $X2 = p1 \cdot \cos(\text{rot2}) \cdot \sin(\text{rot3}) - L \cdot (-\cos(\text{rot3}) \cdot \sin(\text{rot1})) + \cos(\text{rot1}) \cdot \sin(\text{rot2}) \cdot \sin(\text{rot3})) + p2 \cdot (\cos(\text{rot1}) \cdot \cos(\text{rot3}) + \sin(\text{rot1}) \cdot \sin(\text{rot2}) \cdot \sin(\text{rot3}))$ $X3 = -(L \cdot \cos(\text{rot1}) \cdot \cos(\text{rot2})) + p2 \cdot \cos(\text{rot2}) \cdot \sin(\text{rot1}) - p1 \cdot \sin(\text{rot2})$ hence $\tan(\text{Chi}) = X2 / X1$

Parameters

d1: pixel coordinate along the 1st dimation (C convention)

(*type=float or array of them*)

d2: pixel coordinate along the 2nd dimation (C convention)

(*type=float or array of them*)

Return Value

chi, the azimuthal angle in rad

chiArray(*self*, *shape*)

Generate an array of the given shape with chi(i,j) (azimuthal angle) for all elements.

cornerArray(*self*, *shape*)

Generate a 3D array of the given shape with (i,j) (azimuthal angle) for all elements.

cornerQArray(*self*, *shape*)

Generate a 3D array of the given shape with (i,j) (azimuthal angle) for all elements.

delta2Theta(*self*, *shape*)

Generate a 3D array of the given shape with (i,j) with the max distance between the center and any corner in 2 theta

deltaChi(*self*, *shape*)

Generate a 3D array of the given shape with (i,j) with the max distance between the center and any corner in chi-angle

deltaQ(*self*, *shape*)

Generate a 3D array of the given shape with (i,j) with the max distance between the center and any corner in q_vector

diffSolidAngle(*self*, *d1*, *d2*)

calculate the solid angle of the current pixels

solidAngleArray(*self*, *shape*)

Generate an array of the given shape with the solid angle of the current element two-theta(i,j) for all elements.

save(*self*, *filename*)

Save the refined parameters.

Parameters

filename: name of the file where to save the parameters
(*type=string*)

write(*self*, *filename*)

Save the refined parameters.

Parameters

filename: name of the file where to save the parameters
(*type=string*)

sload(*cls*, *filename*)

A static method combining the constructor and the loader from a

Parameters

filename: name of the file to load
(*type=string*)

Return Value

instance of Geometry of AzimuthalIntegrator set-up with the parameter from the file.

load(*self*, *filename*)

Load the refined parameters from a file.

Parameters

filename: name of the file to load
(*type=string*)

read(*self*, *filename*)

Load the refined parameters from a file.

Parameters

filename: name of the file to load
(*type=string*)

getPyFAI(*self*)

return the parameter set from the PyFAI geometry as a dictionary

setPyFAI(*self*, ****kwargs**)

set the geometry from a pyFAI-like dict

getFit2D(*self*)

return a dict with parameters compatible with fit2D geometry

```
setFit2D(self, directDist, centerX, centerY, tilt=0.0, tiltPlanRotation=0.0,
pixelX=None, pixelY=None, splineFile=None)
```

Set the Fit2D-like parameter set: For geometry description see HPR 1996 (14) pp-240

Parameters

direct: direct distance from sample to detector along the incident beam (in millimeter as in fit2d)

tilt: tilt in degrees

tiltPlanRotation: Rotation (in degrees) of the tilt plan around the Z-detector axis * 0deg -> Y does not move, +X goes to Z<0 * 90deg -> X does not move, +Y goes to Z<0 * 180deg -> Y does not move, +X goes to Z>0 * 270deg -> X does not move, +Y goes to Z>0

pixelX, pixelY: as in fit2d they are given in micron, not in meter

centerX, centerY: pixel position of the beam center

splineFile: name of the file containing the spline

```
setChiDiscAtZero(self)
```

Set the position of the discontinuity of the chi axis between 0 and 2pi. By default it is between pi and -pi

```
setChiDiscAtPi(self)
```

Set the position of the discontinuity of the chi axis between -pi and +pi. This is the default behaviour

```
setOversampling(self, iOversampling)
```

set the oversampling factor

```
oversampleArray(self, myarray)
```

```
polarization(self, shape, factor=0.98)
```

Calculate the polarization correction according to the polarization factor:

Parameters

factor: $(I_h - I_v) / (I_h + I_v)$: varies between 0 (no polarization) and 1 (where division by 0 could occur) @return 2D array with polarization correction array (intensity/polarisation)

reset(*self*)

reset most arrays that are cached: used when a parameter changes.

calcfrom1d(*self*, *tth*, *I*, *shape=None*, *mask=None*, *dim1_unit='2th_deg'*)

Computes a 2D image from a 1D integrated profile

Parameters

tth: 1D array with 2theta in degrees

I: scattering intensity @return 2D image reconstructed

set_dist(*self*, *value*)

get_dist(*self*)

set_poni1(*self*, *value*)

get_poni1(*self*)

set_poni2(*self*, *value*)

get_poni2(*self*)

set_rot1(*self*, *value*)

get_rot1(*self*)

set_rot2(*self*, *value*)

get_rot2(*self*)

set_rot3(*self*, *value*)

get_rot3(*self*)

set_wavelength(*self*, *value*)

get_wavelength(*self*)

get_ttha(*self*)

```
set_ttha(self, value)
```

```
del_ttha(self)
```

```
get_chia(self)
```

```
set_chia(self, value)
```

```
del_chia(self)
```

```
get_dssa(self)
```

```
set_dssa(self, value)
```

```
del_dssa(self)
```

```
get_qa(self)
```

```
set_qa(self, value)
```

```
del_qa(self)
```

```
get_pixel1(self)
```

```
set_pixel1(self, pixel1)
```

```
get_pixel2(self)
```

```
set_pixel2(self, pixel2)
```

```
get_splineFile(self)
```

```
set_splineFile(self, splineFile)
```

```
get_spline(self)
```

```
set_spline(self, spline)
```

```
get_correct_solid_angle_for_spline(self)
```

<code>set_correct_solid_angle_for_spline(self, value)</code>
--

Inherited from object

`__delattr__()`, `__format__()`, `__getattribute__()`, `__hash__()`, `__new__()`, `__reduce__()`, `__reduce_ex__()`,
`__setattr__()`, `__sizeof__()`, `__str__()`, `__subclasshook__()`

6.2.2 Properties

Name	Description
dist	
poni1	
poni2	
rot1	
rot2	
rot3	
wavelength	
ttha	2theta array in cache
chia	chi array in cache
dssa	solid angle array in cache
qa	Q array in cache
pixel1	
pixel2	
splineFile	
spline	
correct_SA_spline	
<i>Inherited from object</i>	
<code>__class__</code>	

7 Module pyFAI.geometryRefinement

Date: 23/12/2011

Author: Jerome Kieffer

Contact: Jerome.Kieffer@ESRF.eu

Copyright: European Synchrotron Radiation Facility, Grenoble, France

License: GPLv3+

7.1 Variables

Name	Description
<code>__status__</code>	Value: 'development'
<code>logger</code>	Value: <code>logging.getLogger("pyFAI.geometryRefinement")</code>
<code>ROCA</code>	Value: '/opt/saxs/roca'
<code>__package__</code>	Value: 'pyFAI'

7.2 Class GeometryRefinement

object └

pyFAI.geometry.Geometry └

pyFAI.azimuthalIntegrator.AzimuthalIntegrator └
pyFAI.geometryRefinement.GeometryRefinement

7.2.1 Methods

```
__init__(self, data, dist=1, poni1=None, poni2=None, rot1=0, rot2=0,
rot3=0, pixel1=1, pixel2=1, splineFile=None)
```

`x.__init__(...)` initializes x; see `help(type(x))` for signature

Parameters

`data`: ndarray float64 shape = n, 3 col0: pos in dim0 (in pixels)
col1: pos in dim1 (in pixels) col2: associated tth value (in rad)

Overrides: `object.__init__`

```
residu1(self, param, d1, d2, tthRef)
```

```
residu2(self, param, d1, d2, tthRef)
```

```
refine1(self)
```

```
refine2(self, maxiter=1000000)
```

```
simplex(self, maxiter=1000000)
```

```
anneal(self, maxiter=1000000)
```

```
chi2(self, param=None)
```

```
roca(self)
```

run roca to optimise the parameter set

```
set_dist_max(self, value)
```

```
get_dist_max(self)
```

```
set_dist_min(self, value)
```

```
get_dist_min(self)
```

```
set_poni1_min(self, value)
```

```
get_poni1_min(self)
```

```
set_poni1_max(self, value)
```

```
get_poni1_max(self)
```

```
set_poni2_min(self, value)
```

```
get_poni2_min(self)
```

```
set_poni2_max(self, value)
```

`get_poni2_max(self)``set_rot1_min(self, value)``get_rot1_min(self)``set_rot1_max(self, value)``get_rot1_max(self)``set_rot2_min(self, value)``get_rot2_min(self)``set_rot2_max(self, value)``get_rot2_max(self)``set_rot3_min(self, value)``get_rot3_min(self)``set_rot3_max(self, value)``get_rot3_max(self)`

Inherited from pyFAI.azimuthalIntegrator.AzimuthalIntegrator(Section 3.2)

`makeHeaders()`, `makeMask()`, `reset()`, `save1D()`, `save2D()`, `saxs()`, `setup_LUT()`, `xrpd()`, `xrpd2()`, `xrpd2_histogram()`, `xrpd2_numpy()`, `xrpd2_splitBBox()`, `xrpd2_splitPixel()`, `xrpd_LUT()`, `xrpd_LUT_OCL()`, `xrpd_OpenCL()`, `xrpd_cython()`, `xrpd_numpy()`, `xrpd_splitBBox()`, `xrpd_splitPixel()`

Inherited from pyFAI.geometry.Geometry(Section 6.2)

`__repr__()`, `calcfrom1d()`, `chi()`, `chiArray()`, `chi_corner()`, `cornerArray()`, `cornerQArray()`, `del_chia()`, `del_dssa()`, `del_qa()`, `del_ttha()`, `delta2Theta()`, `deltaChi()`, `deltaQ()`, `diffSolidAngle()`, `getFit2D()`, `getPyFAI()`, `get_chia()`, `get_correct_solid_angle_for_spline()`, `get_dist()`, `get_dssa()`, `get_pixel1()`, `get_pixel2()`, `get_poni1()`, `get_poni2()`, `get_qa()`, `get_rot1()`, `get_rot2()`, `get_rot3()`, `get_spline()`, `get_splineFile()`, `get_ttha()`, `get_wavelength()`, `load()`, `oversampleArray()`, `polarization()`, `qArray()`, `qCornerFunct()`, `qFunction()`, `read()`, `save()`, `setChiDiscAtPi()`, `setChiDiscAtZero()`, `setFit2D()`, `setOversampling()`,

setPyFAI(), set_chia(), set_correct_solid_angle_for_spline(), set_dist(), set_dssa(), set_pixel1(), set_pixel2(), set_poni1(), set_poni2(), set_qa(), set_rot1(), set_rot2(), set_rot3(), set_spline(), set_splineFile(), set_ttha(), set_wavelength(), sload(), solidAngleArray(), tth(), tth_corner(), twoThetaArray(), write()

Inherited from object

__delattr__(), __format__(), __getattr__(), __hash__(), __new__(), __reduce__(), __reduce_ex__(), __setattr__(), __sizeof__(), __str__(), __subclasshook__()

7.2.2 Properties

Name	Description
dist_max	
dist_min	
poni1_min	
poni1_max	
poni2_min	
poni2_max	
rot1_min	
rot1_max	
rot2_min	
rot2_max	
rot3_min	
rot3_max	
<i>Inherited from pyFAI.geometry.Geometry (Section 6.2)</i>	
chia, correct_SA_spline, dist, dssa, pixel1, pixel2, poni1, poni2, qa, rot1, rot2, rot3, spline, splineFile, ttha, wavelength	
<i>Inherited from object</i>	
__class__	

8 Module pyFAI.histogram

Date: 20120916

Author: Jerome Kieffer

8.1 Variables

Name	Description
<code>--package--</code>	Value: 'pyFAI'
<code>--test--</code>	Value: {}

9 Module pyFAI.ocLazim

Date: 03/07/2012

Author: Jerome Kieffer

Contact: jerome.kieffer@esrf.fr

Copyright: 2012, ESRF, Grenoble

License: GPLv3

9.1 Variables

Name	Description
<code>--package--</code>	Value: 'pyFAI'
<code>--test--</code>	Value: {}
<code>lock</code>	Value: <threading.Semaphore object at 0x206f290>
<code>ocl</code>	Value: OpenCL devic...

10 Module `pyFAI.ocl_azim_lut`

Date: 18/10/2012

Author: Jerome Kieffer

Contact: jerome.kieffer@esrf.fr

Copyright: 2012, ESRF, Grenoble

License: GPLv3

10.1 Variables

Name	Description
<code>ocl</code>	Value: OpenCL devic...
<code>__package__</code>	Value: 'pyFAI'

10.2 Class `OCL_LUT_Integrator`

object └─ `pyFAI.ocl_azim_lut.OCL_LUT_Integrator`

10.2.1 Methods

```
__init__(self, lut, devicetype='all', platformid=None, deviceid=None,
checksum=None)
```

`x.__init__(...)` initializes x; see `help(type(x))` for signature

Parameters

lut: array of uint32 - float32 with shape (nbins, lut_size) with indexes and coefficients

checksum: pre - calculated checksum to prevent re - calculating it :)

Overrides: `object.__init__`

```
__del__(self)
```

Destructor: release all buffers

get_nr_threads (<i>self</i> , <i>size</i> =None, <i>ws</i> =None)

calculate the number of threads, multiple of workgroup-size and greater than bins

integrate (<i>self</i> , <i>data</i> , <i>dummy</i> =None, <i>delta_dummy</i> =None, <i>dark</i> =None, <i>flat</i> =None, <i>solidAngle</i> =None, <i>polarization</i> =None)
--

Inherited from object

`__delattr__()`, `__format__()`, `__getattr__()`, `__hash__()`, `__new__()`, `__reduce__()`, `__reduce_ex__()`, `__repr__()`, `__setattr__()`, `__sizeof__()`, `__str__()`, `__subclasshook__()`

10.2.2 Properties

Name	Description
<i>Inherited from object</i>	
<code>__class__</code>	

11 Module pyFAI.peakPicker

Date: 23/12/2011

Author: J\xc3\xa9r\xc3\xebme Kieffer

Contact: Jerome.Kieffer@ESRF.eu

Copyright: European Synchrotron Radiation Facility, Grenoble, France

License: GPLv3+

11.1 Variables

Name	Description
<code>__status__</code>	Value: 'development'
<code>logger</code>	Value: <code>logging.getLogger("pyFAI.peakPicker")</code>
<code>TARGET_SIZE</code>	Value: 1024
<code>__package__</code>	Value: 'pyFAI'

11.2 Class PeakPicker

object └─
 pyFAI.peakPicker.PeakPicker

11.2.1 Methods

<code>__init__(self, strFilename, reconst=False)</code>
<code>x.__init__(...)</code> initializes x; see <code>help(type(x))</code> for signature
Parameters
reconst: shall negative values be reconstucted (wipe out problems with pilatus gaps)
Overrides: <code>object.__init__</code>
<code>gui(self, log=False)</code>
Parameters
log: show z in log scale

load (<i>self</i> , <i>filename</i>)

load a filename and plot data on the screen (if GUI)
--

display_points (<i>self</i>)

onclick (<i>self</i> , <i>event</i>)

readFloatFromKeyboard (<i>self</i> , <i>text</i> , <i>dictVar</i>)

Read float from the keyboard

Parameters

text: string to be displayed

dictVar: dict of this type: {1: [set_dist_min],3: [set_dist_min, set_dist_guess, set_dist_max]}
--

finish (<i>self</i> , <i>filename</i> =None)
--

Ask the 2theta values for the given points
--

contour (<i>self</i> , <i>data</i>)
--

Parameters

data:

massif_contour (<i>self</i> , <i>data</i>)

Parameters

data:

closeGUI (<i>self</i>)

Inherited from object

`__delattr__()`, `__format__()`, `__getattr__()`, `__hash__()`, `__new__()`, `__reduce__()`, `__reduce_ex__()`, `__repr__()`, `__setattr__()`, `__sizeof__()`, `__str__()`, `__subclasshook__()`

11.2.2 Properties

Name	Description
<i>Inherited from object</i>	
<code>__class__</code>	

11.3 Class *ControlPoints*

object —
pyFAI.peakPicker.ControlPoints

This class contains a set of control points with (optionaly) their diffrection 2Theta angle

11.3.1 Methods

`__init__(self, filename=None)`

`x.__init__(...)` initializes `x`; see `help(type(x))` for signature

Overrides: `object.__init__` `exitit`(inherited documentation)

`__repr__(self)`

`repr(x)`

Overrides: `object.__repr__` `exitit`(inherited documentation)

`__len__(self)`

`check(self)`

check internal consistency of the class

`reset(self)`

remove all stored values and resets them to default

`append(self, points, angle=None)`

Parameters

point: list of points

angle: 2-theta angle in radians

`append_2theta_deg(self, points, angle=None)`

Parameters

point: list of points

angle: 2-theta angle in degrees

pop(*self*, *idx=None*)

Remove the set of points at given index (by default the last)

Parameters

idx: position of the point to remove

save(*self*, *filename*)

Save a set of control points to a file

Parameters

filename: name of the file

Return Value

None

load(*self*, *filename*)

load all control points from a file

getList(*self*)

Retrieve the list of control points suitable for geometry refinement

readAngleFromKeyboard(*self*)

Ask the 2theta values for the given points

setWavelength(*self*, *value=None*)

getWavelength(*self*)

Inherited from object

`__delattr__()`, `__format__()`, `__getattribute__()`, `__hash__()`, `__new__()`, `__reduce__()`, `__reduce_ex__()`, `__setattr__()`, `__sizeof__()`, `__str__()`, `__subclasshook__()`

11.3.2 Properties

Name	Description
wavelength	
<i>Inherited from object</i>	
<code>__class__</code>	

11.4 Class Massif

object —
pyFAI.peakPicker.Massif

A massif is defined as an area around a peak, it is used to find neighbouring peaks

11.4.1 Methods

__init__(*self*, *data*=None)

x.__init__(...) initializes *x*; see help(type(*x*)) for signature

Overrides: object.__init__

nearest_peak(*self*, *x*)

@returns the coordinates of the nearest peak

calculate_massif(*self*, *x*)

defines a map of the massif around *x* and returns the mask

find_peaks(*self*, *x*, *nmax*=200, *annotate*=None, *massif_contour*=None, *stdout*=sys.stdout)

All in one function that finds a maximum from the given seed (*x*) then calculates the region extension and extract position of the neighboring peaks.

Parameters

x: seed for the calculation, input coordinates
nmax: maximum number of peak per region
annotate: call back method taking number of points + coordinate as input.
massif_contour: callback to show the contour of a massif with the given index.
stdout: this is the file where output is written by default.

Return Value

list of peaks

initValleySize(*self*)

getValleySize(*self*)

<code>setValleySize(<i>self</i>, <i>size</i>)</code>
--

<code>delValleySize(<i>self</i>)</code>

<code>getBinnedData(<i>self</i>)</code>

@return binned data

<code>getMedianData(<i>self</i>)</code>

<code>getBluredData(<i>self</i>)</code>

<code>getLabeledMassif(<i>self</i>, <i>pattern</i>=None)</code>

Inherited from object

`__delattr__()`, `__format__()`, `__getattr__()`, `__hash__()`, `__new__()`, `__reduce__()`, `__reduce_ex__()`,
`__repr__()`, `__setattr__()`, `__sizeof__()`, `__str__()`, `__subclasshook__()`

11.4.2 Properties

Name	Description
valley_size	Defines the minimum distance between two massifs
<i>Inherited from object</i>	
__class__	

12 Module *pyFAI.reconstruct*

12.1 Variables

Name	Description
<code>--package--</code>	Value: <code>'pyFAI'</code>
<code>--test--</code>	Value: <code>{}</code>

13 Module pyFAI.refinement2D

Date: 23/08/2012

Author: J\xc3\xa9r\xc3\xbbme Kieffer

Contact: Jerome.Kieffer@ESRF.eu

Copyright: European Synchrotron Radiation Facility, Grenoble, France

License: GPLv3+

13.1 Variables

Name	Description
<code>__status__</code>	Value: 'beta'
<code>logger</code>	Value: <code>logging.getLogger("pyFAI.refinement2D")</code>
<code>__package__</code>	Value: 'pyFAI'

13.2 Class Refinement2D

object └─ **pyFAI.refinement2D.Refinement2D**

refine the parameters from image itself ...

13.2.1 Methods

<code>__init__(self, img, ai=None)</code> <hr/> <code>x.__init__(...)</code> initializes x; see <code>help(type(x))</code> for signature Overrides: <code>object.__init__</code>

<code>get_shape(self)</code>

reconstruct (<i>self</i> , <i>tth</i> , <i>I</i>)
--

Reconstruct a perfect image according to 2th / I given in input

Parameters

tth : 2 theta array

I : intensity array

diff_tth_X (<i>self</i> , <i>dx</i> =0.1)

diff_tth_tilt (<i>self</i> , <i>dx</i> =0.1)
--

diff_Fit2D (<i>self</i> , <i>axis</i> ='all', <i>dx</i> =0.1)

scan_centerX (<i>self</i> , <i>width</i> =1.0, <i>points</i> =10)

scan_tilt (<i>self</i> , <i>width</i> =1.0, <i>points</i> =10)
--

scan_Fit2D (<i>self</i> , <i>width</i> =1.0, <i>points</i> =10, <i>axis</i> ='tilt', <i>dx</i> =0.1)
--

Inherited from object

`__delattr__()`, `__format__()`, `__getattr__()`, `__hash__()`, `__new__()`, `__reduce__()`, `__reduce_ex__()`,
`__repr__()`, `__setattr__()`, `__sizeof__()`, `__str__()`, `__subclasshook__()`

13.2.2 Properties

Name	Description
shape	
<i>Inherited from object</i>	
__class__	

14 Module pyFAI.relabel

Date: 20120916

Author: Jerome Kieffer

Contact: Jerome.kieffer@esrf.fr

License: GPLv3+

14.1 Variables

Name	Description
<code>__package__</code>	Value: 'pyFAI'
<code>__status__</code>	Value: 'stable'
<code>__test__</code>	Value: {}

15 Module pyFAI.spline

This piece of software aims to manipulate spline files for geometric corrections of the 2D detectors using cubic-spline

Author: Jérôme Kieffer

Contact: Jerome.Kieffer@esrf.eu

Copyright: European Synchrotron Radiation Facility, Grenoble, France

License: GPLv3+

15.1 Variables

Name	Description
<code>__package__</code>	Value: 'pyFAI'

15.2 Class Spline

This class is a python representation of the spline file. Those files represent cubic splines for 2D detector distortions and make heavy use of fitpack (dierckx in netlib) — A Python-C wrapper to FITPACK (by P. Dierckx). FITPACK is a collection of FORTRAN programs for curve and surface fitting with splines and tensor product splines. See <http://www.cs.kuleuven.ac.be/cwis/research> or <http://www.netlib.org/dierckx/index.html>

15.2.1 Methods

<code>__init__(self, filename=None)</code>
--

this is the constructor of the Spline class, for
--

<code>__repr__(self)</code>

zeros(*self*, *xmin*=0.0, *ymin*=0.0, *xmax*=2048.0, *ymax*=2048.0, *pixSize*=None)

defines a spline file with no (zero) displacement.

Parameters

xmin: minimum coordinate in x, usually zero

(*type*=float)

xmax: maximum coordinate in x (+1) usually 2048

(*type*=float)

ymin: minimum coordinate in y, usually zero

(*type*=float)

ymax: maximum coordinate y (+1) usually 2048

(*type*=float)

zeros_like(*self*, *other*)

defines a spline file with no (zero) displacement with the same shape as the other one given.

Parameters

other: another Spline

(*type*=Spline)

read(*self*, *filename*)

read an ascii spline file from file

Parameters

filename: name of the file containing the cubic spline distortion file

(*type*=string)

comparison(*self*, *ref*, *verbose*=False)

Compares the current spline distortion with a reference

Parameters

ref: another spline file

Return Value

True or False depending if the splines are the same or not

spline2array(*self*, *timing=False*)

calculates the displacement matrix using fitpack bisplev(x, y, tck, dx = 0, dy = 0)

Evaluate a bivariate B-spline and its derivatives. Return a rank-2 array of spline function values (or spline derivative values) at points given by the cross-product of the rank-1 arrays x and y. In special cases, return an array or just a float if either x or y or both are floats.

splineFuncX(*self*, *x*, *y*)

calculates the displacement matrix using fitpack for the X direction

Parameters

x: numpy array representing the points in the x direction

y: numpy array representing the points in the y direction

Return Value

displacement matrix for the X direction

(*type=numpy arrays*)

splineFuncY(*self*, *x*, *y*)

calculates the displacement matrix using fitpack for the Y direction

Parameters

x: numpy array representing the points in the x direction

y: numpy array representing the points in the y direction

Return Value

displacement matrix for the Y direction

(*type=numpy array*)

array2spline(*self*, *smoothing=1000*, *timing=False*)

calculates the spline coefficients from the displacements matrix using fitpack

writeEDF(*self*, *basename*)

save the distortion matrices into a couple of files called basename-x.edf and basename-y.edf

write(*self*, *filename*)

save the cubic spline in an ascii file usable with Fit2D or SPD

Parameters

filename: name of the file containing the cubic spline distortion file
(type=string)

tilt(*self*, *center*=(0.0, 0.0), *tiltAngle*=0.0, *tiltPlanRot*=0.0, *distanceSampleDetector*=1.0, *timing*=False)

The tilt method apply a virtual tilt on the detector, the point of tilt is given by the center

Parameters

center: position of the point of tilt, this point will not be moved.
(type=2tuple of floats)

tiltAngle: the value of the tilt in degrees
(type=float in the range [-90:+90] degrees)

tiltPlanRot: the rotation of the tilt plan with the Ox axis (0 deg for y axis invariant, 90 deg for x axis invariant)
(type=Float in the range [-180:180])

distanceSampleDetector: the distance from sample to detector in meter (along the beam, so distance from sample to center)
(type=float)

Return Value

tilted Spline instance
(type=Spline)

setPixelSize(*self*, *pixelSize*)

sets the size of the pixel from a 2-tuple of floats expressed in meters.

Parameters

pixelSize: *(type=2-tuple of float)*

getPixelSize(*self*)**Return Value**

the size of the pixel from a 2D detector

*(type=2-tuple of floats expressed in meter.)***bin**(*self*, *binning*=None)

16 Module `pyFAI.splitBBox`

16.1 Variables

Name	Description
<code>--package--</code>	Value: <code>'pyFAI'</code>
<code>--test--</code>	Value: <code>{}</code>

17 Module `pyFAI.splitBBoxLUT`

17.1 Variables

Name	Description
<code>--package--</code>	Value: <code>'pyFAI'</code>
<code>--test--</code>	Value: <code>{}</code>

18 Module pyFAI.splitPixel

18.1 Variables

Name	Description
<code>--package--</code>	Value: 'pyFAI'
<code>--test--</code>	Value: {}

19 Module pyFAI.utils

19.1 Functions

timeit(*func*)

gaussian_filter(*input*, *sigma*, *mode*='reflect', *cval*=0.0)

2-dimensional Gaussian filter implemented with FFTw

Parameters

input: input array to filter
(*type=array-like*)

sigma: standard deviation for Gaussian kernel. The standard deviations of the Gaussian filter are given for each axis as a sequence, or as a single number, in which case it is equal for all axes.
(*type=scalar or sequence of scalars*)

mode: {'reflect', 'constant', 'nearest', 'mirror', 'wrap'}, optional The “mode” parameter determines how the array borders are handled, where “cval” is the value when mode is equal to 'constant'. Default is 'reflect'

cval: scalar, optional Value to fill past edges of input if “mode” is 'constant'. Default is 0.0

expand(*input*, *sigma*, *mode*='constant', *cval*=0.0)

Expand array a with its reflection on boundaries

Parameters

a: 2D array

sigma: float or 2-tuple of floats

mode: "constant", "nearest" or "reflect"

cval: filling value used for constant, 0.0 by default

relabel(*label, data, blurred, max_size=None*)

Relabel limits the number of region in the label array. They are ranked relatively to their max(I0)-max(blur(I0))

Parameters

label: a label array coming out of
scipy.ndimage.measurement.label

data: an array containing the raw data

blurred: an array containing the blurred data

max_size: the max number of label wanted @return array like label

averageImages(*listImages, output=None, threshold=0.1, minimum=None, maximum=None, darks=None, flats=None*)

Takes a list of filenames and create an average frame discarding all saturated pixels.

Parameters

listImages: list of string representing the filenames

output: name of the optional output file

threshold: what is the upper limit? all pixel > max*(1-threshold) are discarded.

minimum: minimum valid value or True

maximum: maximum valid value

darks: list of dark current images for subtraction

flats: list of flat field images for division

boundingBox(*img*)

Tries to guess the bounding box around a valid massif

Parameters

img: 2D array like

Return Value

4-tuple (d0_min, d1_min, d0_max, d1_max)

removeSaturatedPixel(*ds*, *threshold*=0.1, *minimum*=None, *maximum*=None)

Parameters

ds: a dataset as ndarray
threshold: what is the upper limit? all pixel > max*(1-threshold) are discarded.
minimum: minimum valid value (or True for auto-guess)
maximum: maximum valid value

Return Value

another dataset

binning(*inputArray*, *binsize*)

Parameters

inputArray: input ndarray
binsize: int or 2-tuple representing the size of the binning

Return Value

binned input ndarray

unBinning(*binnedArray*, *binsize*)

Parameters

binnedArray: input ndarray
binsize: 2-tuple representing the size of the binning

Return Value

unBinned input ndarray

19.2 Variables

Name	Description
logger	Value: logging.getLogger("pyFAI.utils")
timelog	Value: logging.getLogger("pyFAI.timeit")
__package__	Value: 'pyFAI'

Index

- pyFAI (*package*), 2
 - pyFAI._geometry (*module*), 3
 - pyFAI.azimuthalIntegrator (*module*), 4–20
 - pyFAI.azimuthalIntegrator.AzimuthalIntegrator (*class*), 4–20
 - pyFAI.bilinear (*module*), 21
 - pyFAI.detectors (*module*), 22–30
 - pyFAI.detectors.Detector (*class*), 22–24
 - pyFAI.detectors.Fairchild (*class*), 28–29
 - pyFAI.detectors.FReLoN (*class*), 29–30
 - pyFAI.detectors.Pilatus (*class*), 24–25
 - pyFAI.detectors.Pilatus1M (*class*), 25–26
 - pyFAI.detectors.Pilatus2M (*class*), 26–27
 - pyFAI.detectors.Pilatus6M (*class*), 27–28
 - pyFAI.geometry (*module*), 31–43
 - pyFAI.geometry.Geometry (*class*), 31–43
 - pyFAI.geometryRefinement (*module*), 44–47
 - pyFAI.geometryRefinement.GeometryRefinement (*class*), 44–47
 - pyFAI.histogram (*module*), 48
 - pyFAI.ocl_azim (*module*), 49
 - pyFAI.ocl_azim_lut (*module*), 50–51
 - pyFAI.ocl_azim_lut.OCL_LUT_Integrator (*class*), 50–51
 - pyFAI.peakPicker (*module*), 52–57
 - pyFAI.peakPicker.ControlPoints (*class*), 53–55
 - pyFAI.peakPicker.Massif (*class*), 55–57
 - pyFAI.peakPicker.PeakPicker (*class*), 52–53
 - pyFAI.reconstruct (*module*), 58
 - pyFAI.refinement2D (*module*), 59–60
 - pyFAI.refinement2D.Refinement2D (*class*), 59–60
 - pyFAI.relabel (*module*), 61
 - pyFAI.spline (*module*), 62–66
 - pyFAI.spline.Spline (*class*), 62–66
 - pyFAI.splitBBox (*module*), 67
 - pyFAI.splitBBoxLUT (*module*), 68
 - pyFAI.splitPixel (*module*), 69
 - pyFAI.utils (*module*), 70–72
 - pyFAI.utils.averageImages (*function*), 71
 - pyFAI.utils.binning (*function*), 72
 - pyFAI.utils.boundingBox (*function*), 71
 - pyFAI.utils.expand (*function*), 70
 - pyFAI.utils.gaussian_filter (*function*), 70
 - pyFAI.utils.relabel (*function*), 70
 - pyFAI.utils.removeSaturatedPixel (*function*), 71
 - pyFAI.utils.timeit (*function*), 70
 - pyFAI.utils.unBinning (*function*), 72